

Important information regarding iPad cases

It is important that you order your keyguard to fit the case you intend to use it with. On the following pages you will find information about the cases that we currently cut the keyguards to fit. Please be sure to review the information on your respective keyguard before you order. We will make your keyguard according to your instructions, and for the case you order. Once the keyguard has been cut it may not be possible to change it, with very few exceptions. If you fail to order the keyguard to fit the correct case it is very likely that it will not fit correctly and you will need to order the correct keyguard.

Please keep in mind that the keyguards are not designed to attach to any particular case though some cases do offer attachment options that are available for an additional charge. For cases without an attachment option there are other possibilities to keep the keyguard in place such as suction cups and/or elastic straps. Both of these are available on our website for an additional charge.

Ultimate II by RJ Cooper

Suction cups will not work for keyguards made to fit this case, with a few exceptions. You may use the optional Strap Attachment (additional purchase), item #SA-337 with this case.



Otterbox Defender

The Otterbox Defender case for the iPad 1 is similar to the Defender case for the iPad 2/3/4, shown here. Keyguards for this case will not support suction cups. Strap Attachments, item #SA-337(available as an optional, additional purchase), is available for most keyguards made for this size case.




The Otterbox Defender for the iPad 2, 1st generation (shown below), will support the attachment of suction cups (strap attachment may also be used with this case). A good indicator of which case you have is whether or not the camera lens is surrounded by the case, as shown by the arrow above, or open, as shown by the arrow below. If open, as below, the case is a 1st generation and the keyguard will be cut similar to the image on the right.



Images are not to scale

AMDi iAdapter

The AMDi iAdapter case is available for all versions of iPads. These keyguards do not allow for the attachment of suction cups. There are several options to attach the keyguard. One, you can use the strap attachment, item #SA-337. Two, you can use the case attachment option, item #ATT-1114 or Three, you may use the AMDi iAdapter Attachment, item #AMDi3. These are all options and require separate purchase from the keyguard.

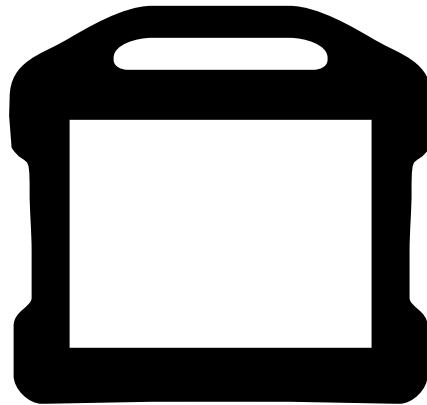

The Case Attachment option, item #ATT-1114 is made by cutting the keyguard and then cutting a second piece of acrylic slightly larger than the keyguard. The two are glued together to make one piece which, when viewed on it's edge, will look like the image below. The left and right edge will rest on the iAdapter case with the keyguard resting on the iPad. You will need to use Velcro, not included, to keep the keyguard in place.



The AMDi iAdapter Attachment, item #AMDi3, shown below, is cut from one piece of acrylic with a portion that has been cut away. The “ledges” on either side will slide in between the case and the iPad. You will need to slightly loosen the screws on the back of the case then slip one side into the gap between the iPad and the iAdapter case. You will then GENTLY flex the keyguard and slide the other side into the gap. You may then tighten the screws.



Images are not to scale

Big Grips

Keyguards for the Big Grips case will allow for the use of suction cups. You may also use the optional Strap Attachment (additional purchase), item #SA-337 with this case.



Griffin Survivor


Keyguards for the Griffin Survivor will not allow for the use of suction cups. You may use the optional Strap Attachment (additional purchase), item #SA-337 with this case. NOTE: Some keyguards, such as ones made for LAMP Words for Life, Speak for Yourself, etc., will not allow for straps due to the number of openings required by the app.



Images are not to scale

Gumdrop Case

Keyguards made for the Gumdrop case will not allow for the use of suction cups. You may use either the Strap Attachment, item #SA-337 or the Gumdrop Attachment (shown below), item #GA-147. Both of these are optional and are an additional charge. The Strap Attachment may be added at any time, the Gumdrop Attachment must be ordered at the same time as your keyguard, it cannot be added later. The Gumdrop Attachment will require removing the screen protector insert from the case. There are four tabs, two on either side of the keyguard, that insert where the screen protector was. You will insert one side then GENTLY flex the keyguard and insert the other side

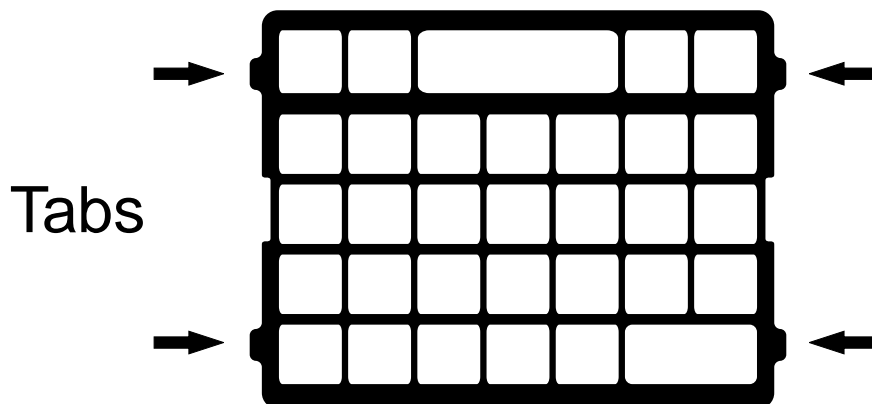

Images are not to scale

XGR360

The Xgear 360 will allow for the use of suction cups. You may also use the optional Strap Attachment(additonal purchase), item #SA-337 with this case.



GoNow Case

The GoNow case will not allow for the use of suction cups. You may use the optional Strap Attachment(additonal purchase), item #SA-337 with this case.



Images are not to scale

KaysCase Kidbox Monoprice Kidz Cover i-Blason Armor Box

Suction cups will not work for keyguards made to fit this case, with a few exceptions. You may use the optional Strap Attachment(additonal purchase), item #SA-337 with this case.



M-Edge Super Shell

Keyguards for this case can use both suction cups, item #CUP, and also straps, item #SA-337. Both items are optional and available for an additional charge.



Lifeproof Nuud

Suction cups will not work for keyguards made to fit this case, with a few exceptions. You may use the optional Strap Attachment(additonal purchase), item #SA-337 with this case.



Grip Case

Keyguards for this case can use both suction cups, item #CUP, and also straps, item #SA-337. Both items are optional and available for an additional charge.

